



D-LUX+ CORDLESS LED CURING LIGHT NEW

D-LUX+ is a high performance dual-wavelength Cordless LED Curing Light. Built and designed with the highest quality and durability, its dual-wavelength can efficiently, effectively, and uniformly polymerize all light-cured dental materials. At its maximum light intensity output of 2,400mW/cm², it can complete all restorations safely and efficiently.




Features


- A high intensity output of **2,400mW/cm²** and 10W LED offers short curing time and efficient polymerization
- 360° swivel light probe and silicone ring tip ensure patient comfort
- **Dual-wavelength** with 4 LED lamps can efficiently polymerize all dental materials
- Broadband spectrum of **385nm~515nm**
- Five (5) versatile curing modes allow flexible curing time and power intensity including 1second curing and Orthodontic curing mode
- **Wireless battery charging system** eliminates any potential charging contact problems
- High capacity 2,600mAh Li-ion battery offers longer battery life
- Built-in Light Intensity Measurement Sensor provides peace of mind so clinicians can always verify if the intensity is within normal range
- Designed ergonomically, the handpiece is light and elegant so it can easily reach all curing sites and offers maximum comfort for both patient and clinician

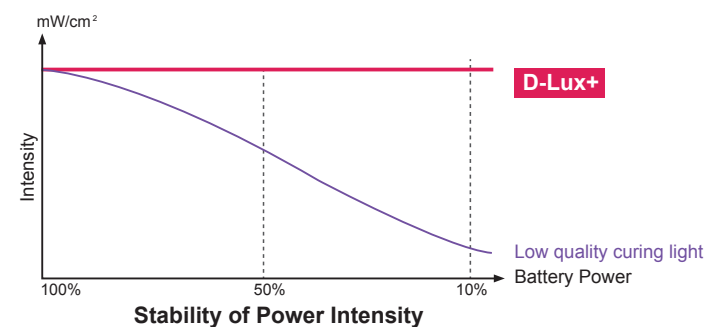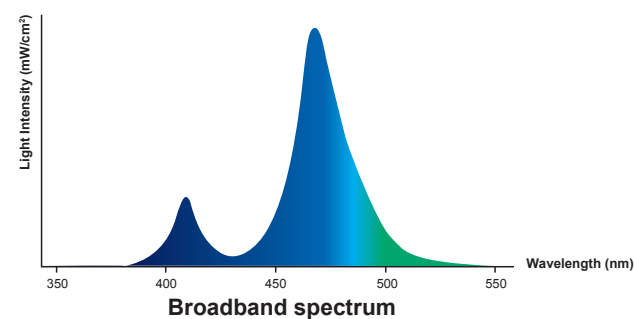


5 Curing Modes

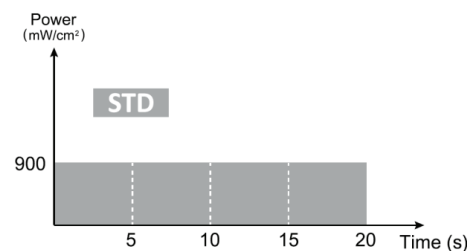


STD (Standard):
900mW/cm²(±10%) light intensity for 5, 10, 15, or 20 seconds

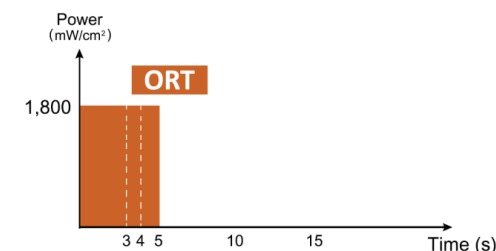


ORT (Orthodontic):
1,800mW/cm²(±10%) light intensity for 3, 4, 5 seconds or P (pulse)

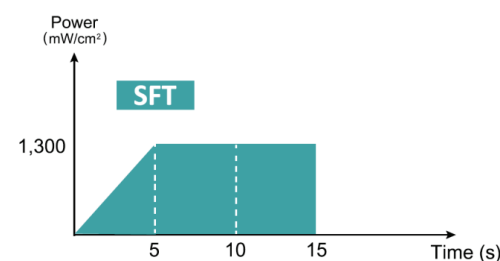


SFT (Soft start):
Light intensity gradually increases to 1,300mW/cm²(±10%) for 10 or 15 seconds

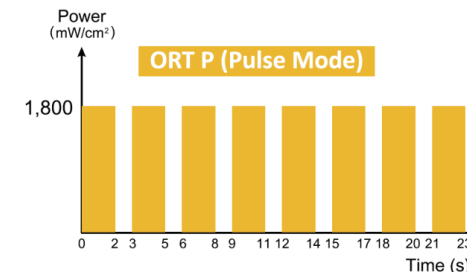

ORT P (Orthodontic Pulse Mode):
1,800 mW/cm² (±10%) of light turns on for 2 seconds and then off for 1 second eight times

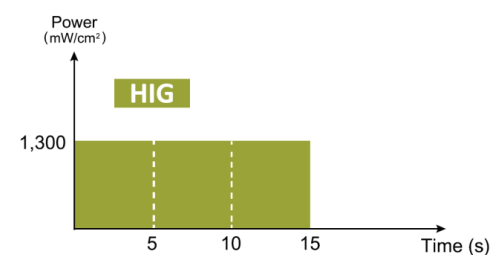

HIG (High Power):
1,300mW/cm²(±10%) light intensity for 5, 10, or 15 seconds

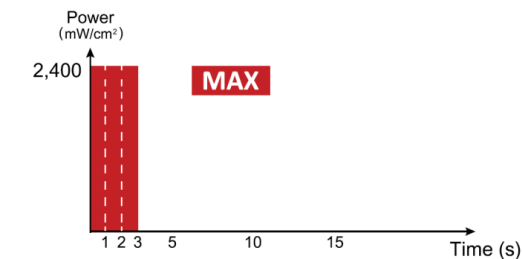

MAX (Max Power):
2,400mW/cm²(±10%) light intensity for 1, 2, or 3 seconds

Product Features

Wireless battery charging system

Dual-wavelength with 4 LED lamps, Broadband Spectrum of 385nm~515nm

Built-in light intensity sensor

Package

- ITEM #4008-1001 Regular Kit
 - D-Lux+ Handpiece & Charger
 - Light Probe
 - Light Protector
 - Disposable Sheaths(200ea/Box)
 - Power Supply(Adapter), Power Cord

Accessories

- ITEM #4008-1110 Handpiece
- ITEM #4008-1120 Charger
- ITEM #4008-1130 Light Probe
- ITEM #4008-1140 Battery Pack
- ITEM #4008-1150 Light Protector
- ITEM #4008-1160 Power Supply(Adapter)
- ITEM #4008-1161 Power Cord
- ITEM #4008-1170 Disposable Sheaths(200ea/Box)